


Home

To start word prediction:

- Click on one of the suggested words.
- OR Start the recording. When it finishes, your speaking is transformed to written words and labeled with images

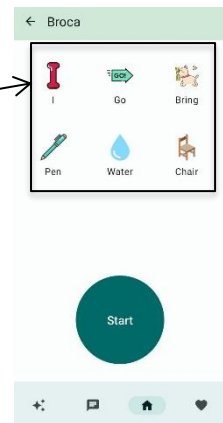

After that:

- A row with suggested words to continue your sentence will appear. Click on the desired word to add it in the sentence. New suggested words will appear automatically.
- Click Start Button to complete the sentence
- Click Reset Button to reset the sentence
- Click on the Arrow to process the sentence



Process sentences

After clicking on the "Process" button, a list of suggested sentences based on your words is generated. The first given sentence is the original non-processed sentence.

- Click on Audio Button to listen the sentence
- Click on Favorite Button to save the sentence and reuse it later


Saved Sentences

This sentence will be available when needed, you will be able to read it or hear it by clicking on the Audio Button.

You can save frequently used sentences by 2 means:

- Click on Favorite Button on a sentence in the Process screen
- Add a sentence directly writing it down using Text Input

To delete Saved Sentence, re-click on its Favorite Button.

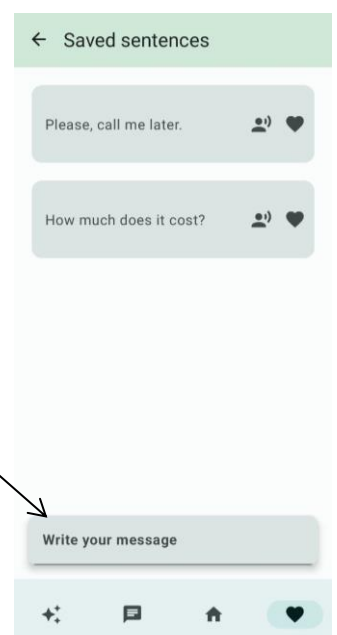

ChatBot



To practice your language skills, you may try to speak with chatbot. The chatbot can also help you to interact with other people. To transform into Speech the chatbot response, click on Audio button at the end of the message.

Practice



Broca has the same functionalities as other Aphasia Apps in terms of therapeutic exercises.

Helps to recreate links between images and words.

2 types of training are available:

- Association of written words to the image
- Association of spoken words to the image



Click on one of the four options to submit the answer, the button will become green if the answer is right, red otherwise.

Click on the Next button to get the next word to guess.

